Instructions for Use

Anti-tTG IgA

IVD

CE

REF

HYE-5894

40 tests

DRG Instruments GmbH, Germany
Frauenbergstraße 18, D-35039 Marburg
Phone: +49 (0)6421-1700 0, Fax: +49 (0)6421-1700 50
Website: www.drg-diagnostics.de
E-mail: drg@drg-diagnostics.de

Distributed by:

DRG International, Inc., USA
841 Mountain Ave., Springfield, NJ 07081
Phone: (973) 564-7555, Fax: (973) 564-7556
Website: www.drg-international.com
E-mail: corp@drg-international.com

ASSAY PROTOCOL BARCODE (APPB)

(Version 2.61 Software or later / ab Softwareversion 2.61 oder höher / Versione Software 2.61 o superiore / Versión de Software 2.61 o posterior / Version de logiciel 2.61 et supérieure / Wersja oprogramowania 2.61 lub późniejsza)

HYE-5894 – v2.61.1

The barcode must be used to install the assay protocol into the DRG:HYBRID-XL software via the SCAN NEW LOT page.
Der Barcode muss in dem Menü „NEUES LOT SCANNEN“ eingelesen werden, um das Protokoll in der DRG:HYBRiD-XL-Software zu installieren.

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Please refer to section 3: **Routine Procedures: “Installing a new assay product”** of the User Manual v2.60 or later.

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Contents / Inhaltsverzeichnis / Contenuti / Contenido / Contenu / Spis treści

1	INTRODUCTION	4
2	PRINCIPLE OF THE TEST	4
3	WARNINGS AND PRECAUTIONS	4
4	REAGENTS	4
5	SPECIMEN COLLECTION AND PREPARATION	5
6	ASSAY PROCEDURE	6
7	QUALITY CONTROL	6
8	EXPECTED NORMAL VALUES	6
9	LIMITATIONS OF USE	6
10	PERFORMANCE CHARACTERISTICS	6
11	METHOD COMPARISON	7
12	LEGAL ASPECTS	7

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1 INTRODUCTION

1.1 Intended Use

The DRG:HYBRiD-XL Anti-tTG IgA is an enzyme immunoassay for the quantitative *in vitro diagnostic* measurement of IgA class antibodies against anti-tissue-transglutaminase (anti-tTG IgA) in serum or EDTA plasma.

Only for use with the DRG:HYBRiD-XL Analyzer.

1.2 Summary and Explanation

Celiac disease (CD; synonym: gluten-sensitive enteropathy or gluten intolerance) is a systemic autoimmune disease characterized by an immunological reaction after uptake of gluten (prolamines and glutelins) from cereal grains (e.g. wheat, oat, barley or rye) (1). The disease occurs in genetically susceptible individuals and is associated with HLA alleles DQ2 or DQ8 (1). The worldwide prevalence of CD is 0.5-1% of the general population (3), but may be 2-5% in some populations (4,16). Furthermore, risk factors like early and massive gluten exposure, bowel infections, or drugs may contribute to CD (3). The disease causes morphological changes of the upper small intestine with atrophy of the villi of the mucous membrane. As a consequence, the resorption surface of the small intestine is reduced, leading to malabsorption, secondary pancreatic insufficiency, or chronic vitamin deficiency (5).

It has been known for many years that elevated levels of gluten-specific antibodies occur in the sera of celiac patients. Furthermore, autoantibodies directed against endomysium of the smooth muscle are specifically associated with this disease (6-9). Among these, tissue transglutaminase (tTG) has been identified as the predominant endomysium CD autoantigen (8). Celiac disease constitutes an indication for further autoimmune diagnostic testing, since 30% of cases suffer from autoimmune diseases like juvenile diabetes (10), autoimmune thyroiditis (Hashimoto's and Graves' disease), and autoimmune hepatitis (10,11,14,17), or various forms of arthritis (12-13).

The only effective therapy for celiac disease is a life-long, strict gluten-free diet because even small amounts of gluten can induce pathological conditions (1,5-6).

Tissue transglutaminase is a calcium-dependent enzyme which exists particularly in the intestine (5). tTG deamidates gluten peptides, especially gliadin, resulting in glutamic acid residues, thereby promoting the binding of gluten peptides to HLA-DQ2 or HLA-DQ8 and potentiating the inflammatory T-cell reaction. The activated T cells in the lamina propria and within the epithelium are cytotoxic, with effects including apoptosis of enterocytes, atrophic remodelling of the mucosa, and malabsorption. Furthermore, antibodies to gliadin peptides and tTG are formed, which either accumulate in the tissue or enter the bloodstream (1).

The serological determination of celiac-specific parameters is performed using immunoassays. The parameters include antibodies to the tissue transglutaminase, deamidated gliadin peptides or endomysium. Anti-tTG IgA-ELISA is recommended as the primary serological screening test for celiac disease. The total IgA level should be measured before the test, because more than 2% of patients with celiac disease have a selective IgA deficiency, so that IgA antibody tests may be negative even in the phase of active disease. If IgA deficiency is found, a test for IgG autoantibodies against deamidated gliadin peptides (IgG anti-DGP) is recommended (1,14-15). By serological examination, a burdensome biopsy can be avoided. The test is also well suited for monitoring and controlling the gluten-free diet, as the antibody level drops on the diet and increases with renewed gluten intake. (1,5).

2 PRINCIPLE OF THE TEST

The DRG:HYBRiD-XL Anti-tTG IgA Kit is a solid phase enzyme-linked immunosorbent assay (ELISA) based on the **indirect immunoassay principle**.

The wells (ACW) of the reagent cartridges are coated with a tissue transglutaminase (tTG).

An aliquot of automatically pre-diluted patient sample containing endogenous anti-tTG IgA antibodies is incubated in the coated well. After incubation the unbound sample components are washed off. Thereafter, enzyme conjugate, which is an anti-human IgA antibody conjugated with horseradish peroxidase, is incubated in the coated well. After incubation the unbound conjugate is washed off.

The amount of bound peroxidase conjugate is proportional to the concentration of anti-tTG IgA in the sample.

Having added the substrate solution, the intensity of colour developed is proportional to the concentration of anti-tTG IgA in the patient sample.

3 WARNINGS AND PRECAUTIONS

- This kit is for *in vitro diagnostic* use only.
For professional use only.

- This kit can only be used in combination with the DRG:HYBRiD-XL Analyzer
- Before starting the assay, read the instructions completely and carefully. Use the valid version of the package insert provided with the kit. Be sure that everything is clear and understood.
- Do not remove, exchange, discard or damage any of the barcode labels provided with each kit and its components. All barcodes build an integral system for the kit lot.**
- Respect the general safety measures for use of laboratory reagents.
- All reagents of this test kit which contain human serum or plasma have been tested and confirmed negative for HIV I/II, HBsAg and HCV by FDA approved procedures. All reagents, however, should be treated as potential biohazards in use and for disposal.
- Never pipet by mouth and avoid contact of reagents and specimens with skin and mucous membranes.
- Do not smoke, eat, drink or apply cosmetics in areas where specimens or kit reagents are handled.
- Wear appropriate disposable gloves when handling specimens and reagents. Microbial contamination of reagents or specimens may cause false results.
- Handling should be done in accordance with the procedures defined by an appropriate national biohazard safety guideline or regulation.
- Do not use reagents beyond expiry date as shown on the kit labels.
- Unused reagent cartridges must be stored at 2 °C to 8 °C in the sealed foil pouch with desiccant provided.
- Optimal test results are only obtained when using calibrated pipettes.
- Do not mix or use components from kits with different lot numbers. It is advised not to interchange reagent cartridges of different kits even of the same lot. The kits may have been shipped or stored under different conditions and the binding characteristics of the wells in the reagent cartridges may differ slightly.
- Some reagents contain Proclin 300, BND and/or MIT as preservatives. In case of contact with eyes or skin, flush immediately with water.
- TMB substrate has an irritant effect on skin and mucosa. In case of possible contact, wash eyes with an abundant amount of water and skin with soap and plenty of water. Wash contaminated objects before reusing them. If inhaled, take the person to open air.
- Chemicals and prepared or used reagents have to be treated as hazardous waste according to the national biohazard safety guidelines or regulations.
- For information on hazardous substances included in the kit please refer to Safety Data Sheets. For professional users the Safety Data Sheet for this product is available upon request directly from DRG.

4 REAGENTS

4.1 Reagents provided

4.1.1 Reagent Cartridges

40 pieces containing the following:

- **Antibody Coated Well (ACW)**
coated with tissue transglutaminase (tTG).
- **Sample Diluent**, 300 µL,
Contains non-mercury preservative.
- **Enzyme Conjugate**, 170 µL,
Anti-human IgA antibody conjugated with horseradish peroxidase;
Contains non-mercury preservative.
- **Substrate Solution**, 270 µL
Tetramethylbenzidine (TMB).

4.1.2 Additional Cartridges (red)

40 pieces containing the following:

- **Sample Diluent**, 300 µL,
Contains non-mercury preservative.

4.1.3 Re-Calibrator 1 & 2

2 vials, 1.5 mL each, ready to use;

For re-calibration of the quantitative DRG:HYBRiD-XL Anti-tTG IgA test.

Concentrations are lot-specific.

Contain non-mercury preservative.

4.1.4 Control 1 & 2

2 vials, 1.5 mL each, ready to use;

For control values and ranges please refer to the bar code on vial label or to the QC-Datasheet.

Contain non-mercury preservative.

4.2 Materials required but not provided

- General needed laboratory equipment

- Ultra-pure water

DRG recommends to use Clinical Laboratory Reagent Water (CLRW) according to CLSI guideline 3C-A4 with the following specifications:

Resistivity at 25 °C [MΩ·cm]: > 10

Conductivity at 25 °C [$\mu\text{S}\cdot\text{cm}^{-1}$]: < 0.1

Total Organic Carbon/p.p.b.[$\mu\text{g}/\text{L}$]: < 50

Colloids [$\mu\text{g}/\text{mL}$]: <0.05

- **REF** HYB-5252 DRG:HYBRID-XL Analyzer

- **REF** HYI-5392: System Solution 5L, 5000 mL;

(Instrument Feed Water according to CLSI guideline 3C-A4 with the following specification can also be used:

Resistivity at 25 °C [MΩ·cm]: > 1

Conductivity at 25 °C [$\mu\text{S}\cdot\text{cm}^{-1}$]: < 1

Total Organic Carbon/p.p.b.[$\mu\text{g}/\text{L}$]: < 200

Colloids [$\mu\text{g}/\text{mL}$]: <0.1)

- **REF** HYI-5394: Wash Buffer, 40x concentrate, 25 mL

- **REF** HYI-5395: Needle Cleaning Solution, 30 mL. Cleaning solution for the pipetting needle (daily and weekly maintenance, see also user manual)

- **REF** HYI-5387: Cuvettes, (2 x 360 pieces)

For use of the *Secondary Sample Holder* for secondary tubes the following tubes are required:

- **REF** HYI-5391: Sample Tubes (Secondary), 2500 pcs.

4.3 Storage Conditions

All kit components should be stored at 2 °C to 8 °C to ensure product performance until the defined expiry date.

When stored at 2 °C to 8 °C, **unopened kits** will retain reactivity until expiration date. Do not use reagents beyond this date.

- Cartridges (stored at 2 °C to 8 °C) in the supplied and unopened zip/foil bags will retain reactivity until expiration date.
- Unopened Re-Calibrators and Controls (stored at 2 °C to 8 °C) will retain reactivity until expiration date.

Opened reagents and the reagent cartridges must be stored at 2 °C to 8 °C.

Once the plastic bag has been opened, care should be taken to tightly close it again along with the supplied desiccant bag.

Immediately after end of each run the Re-Calibrator and Control vials have to be removed from the instrument, tightly capped and stored at 2 °C to 8 °C.

- Unused cartridges in opened zip/foil bags (stored at 2 °C to 8 °C) will retain reactivity until expiration date, if stored as described above.
- Pierced or open cartridges must be disposed of immediately.
- Opened Re-Calibrators and Controls (stored at 2 °C to 8 °C) will retain reactivity for 8 weeks.

4.3.1 On-board Stability

For Re-Calibrators and Controls the on-board stability has been evaluated under controlled laboratory conditions at room temperature (20 °C to 25 °C).

Due to the differences in laboratory environmental conditions and reagent volumes, the on-board stability may deviate from the declared value.

On-board stability:	24 h
---------------------	------

4.4 Reagent Preparation

Bring all reagents, such as controls and re-calibrators, to room temperature (20 °C to 25 °C) prior to use. Reagent Cartridges can be used directly after storage in the refrigerator.

Wash Buffer (not included in the kit)

For Wash Buffer (1x) dilute 25 mL of *Wash Buffer* (40x) with 975 mL ultra-pure water to a final volume of 1000 mL.

The diluted Wash Buffer (1x) is stable for 2 weeks at room temperature.

4.5 Disposal of the Kit

The disposal of the kit and all used materials/reagents must be performed according to the national regulations. Special information for this product is given in the Safety Data Sheet.

4.6 Damaged Test Kits

In case of any damage to the test kit or components, DRG must be informed in writing, at the latest one week after receiving the kit.

Damaged single components should not be used for a test run. They have to be stored until a final solution has been found. After this, they should be disposed of according to the official regulations.

5 SPECIMEN COLLECTION AND PREPARATION

Serum or EDTA plasma can be used in this assay.

A minimum of 72 μL of sample is needed for one determination.

This includes 12 μL sample and 60 μL dead volume.

A minimum of 210 μL of Control and Re-Calibrator (ready-to-use) is needed for one determination. This includes 150 μL Control and Re-Calibrator and 60 μL dead volume.

Attention:

- This test was not verified with blood collection tubes of all available manufacturers.
- Sample Collection Systems of some manufacturers may contain different materials which in isolated cases could affect the test results.
- If primary tubes for sample collection are used, please follow the instructions of the manufacturer.
- Do not use haemolytic, icteric or lipaemic specimens.
- Samples containing precipitates have to be centrifuged prior to the test run.
- Do not use heat inactivated samples.
- Do not use standards or external controls stabilized with azide.

5.1 Specimen Collection**Serum:**

Collect blood by venipuncture (e.g. Sarstedt Monovette for serum), allow to clot, and separate serum by centrifugation at room temperature. Do not centrifuge before complete clotting has occurred. Patients receiving anticoagulant therapy may require increased clotting time.

Plasma:

Whole blood should be collected into centrifuge tubes containing anti-coagulant (e.g. Sarstedt Monovette with the appropriate plasma preparation) and centrifuged immediately after collection.

5.2 Specimen Storage and Preparation

Specimens should be capped and may be stored for up to 7 days at 2 °C to 8 °C prior to performing the assay.

Specimens stored for a longer time (up to 9 months) should be frozen only once at -20 °C prior to the assay. Thawed samples should be inverted several times prior to testing.

5.3 Specimen Dilution**5.3.1 Fully Automated Pre-dilution of Samples**

Each sample will be **pre-diluted fully automated 1:200** using the additional cartridge (red) included in the kit (1 cartridge per 1:200 dilution).

5.3.2 Manual Sample Dilution

If in an initial assay, a specimen is found to contain more than the highest standard, the specimens can be diluted with *Sample Diluent** and measured again as described in Assay Procedure.

For the calculation of the concentrations this dilution factor has to be taken into account.

Example:

a) dilution 1:10: 10 μL sample + 90 μL *Sample Diluent* (mix thoroughly)

* *Sample Diluent for manual dilution is not included in this kit, but can be ordered on request (REF HYE-5894-DIL, 20 mL).*

6 ASSAY PROCEDURE

6.1 General Remarks

- All reagents, such as controls and re-calibrators and specimens, must be allowed to come to room temperature (20 °C to 25 °C) before use.
All reagents and samples must be mixed without foaming.
Reagent Cartridges can be used directly after storage in the refrigerator.
- Samples, controls and re-calibrators should be measured within 2 hours in order to avoid possible evaporation effects.
- The **Secondary Sample Holder** (HY1-5437) for secondary tubes has the capacity for a maximum of 20 samples including controls and re-calibrators. They all have to be pipetted into the secondary tubes, and the respective barcodes of control/re-calibrator vials and, if available, the sample barcodes have to be read with the external barcode scanner.

6.2 Test Procedure

- **The total assay time for DRG:HYBRiD-XL Anti-tTG IgA is 120 minutes.**
- To ensure proper operation of the test, the instructions in the user manual for the DRG:HYBRiD-XL should strictly be followed.
- All test specific information required for the correct operation is included in the respective barcodes of the reagents.
Take care not to damage these bar codes!
- **It is recommended to tap the bottom of the Cartridge Segments containing the reagent cartridges once on the bench before placing them on the rotor. This is to avoid foam and adhering of the liquid on the sealing of the reagent cartridge.**
- Place reagent cartridges on the rotor of the unit. The heating to 37 °C incubation temperature is performed automatically in the unit.

6.3 Calibration

Traceability:

This method was standardized against the DRG t-Transglutaminase-Ab IgA ELISA (EIA-4259).

Each DRG:HYBRiD-XL reagent contains a barcode with the specific information for recalibration of the reagent lot. The Master Curve is printed as a 2-D barcode on the outer label of the kit package and on the QC-Datasheet and has to be scanned with the external barcode scanner prior to the first use of the respective kit lot.

Recalibration is recommended:

- if a new kit lot is used. Each new lot should be verified by running the kit internal re-calibrators and controls before routine use.
- if one or both assay controls are found outside the specified range.
- after 4 weeks of use of the same reagent kit on the unit.

6.4 Calculation of Results

The analyte concentrations are calculated automatically by the DRG:HYBRiD-XL's system software.

7 QUALITY CONTROL

It is recommended to use control samples according to state and federal regulations. The use of control samples is advised to assure the day to day validity of results.

It is also recommended to participate in national or international Quality Assessment programs in order to ensure the accuracy of the results. Apply appropriate statistical methods for analysing control values and trends. If the results of the assay do not agree with the established acceptable ranges of control materials, patient results should be considered invalid. In this case, please check the following: expiration dates and storage conditions of reagents, operational reliability of the analyser. In addition, it is indicated to perform a Recalibration.

In case of further questions please contact your local distributor or DRG directly.

7.1 Internal Controls

For Quality Control it is necessary to use the two internal controls provided with each kit.

Acceptance ranges for both internal controls (*Control 1 & 2*) were established by the manufacturer and are summarized in the QC-Datasheet added to the kit. Note that the expected values and acceptance ranges stated in the QC-Datasheet always refer to the current kit lot.

Internal controls should be run in single determination:

- on a routine basis (e.g. once per 24 h)
- if re-calibration is required (if one or both internal controls are out of range)

- if a new kit lot is used (in order to avoid any negative impact on the kit performance by improper transport or to detect improper storage during transport).

7.2 External Controls

Use controls at both normal and pathological levels.

The control intervals and control ranges for external controls should be adapted to the individual requirements of each laboratory. All results must be within the defined limits.

Each laboratory should establish corrective measures to be taken if values of external controls are not found in the acceptance range.

8 EXPECTED NORMAL VALUES

It is strongly recommended that each laboratory should determine its own normal and pathological values.

In a study conducted with apparently healthy individuals, using the DRG:HYBRiD-XL Anti-tTG IgA the following data were observed:

Population	Males	Females	Total
n	60	59	119
Range (min. - max.) (U/mL)	0.2 - 3.9	0.2 - 2.3	0.2 - 3.9
Mean (U/mL)	1.0	0.7	0.8
2.5 th - 97.5 th Percentile (U/mL)	0.3 - 2.5	0.2 - 2.0	0.3 - 2.1
Median (U/mL)	0.8	0.5	0.7

8.1 Cut-off Value

The Cut-off Value of the DRG:HYBRiD-XL Anti-tTG IgA was determined by measuring the values of 119 apparently healthy individuals and 60 patients with coeliac disease and a positive anti-tTG IgA result (obtained by the reference method ELISA (DRG, EIA-4259)).

	CD patients
n	60
Range (min. - max.) (U/mL)	2.6 - > 100
Mean (U/mL)	54.4
2.5 th - 97.5 th Percentile (U/mL)	2.89 - > 100
Median (U/mL)	43.8

The Cut-off value was defined as 3.1 U/mL.

Using a Cut-off value of 3.1 U/mL, the specificity was calculated to 99.2% and sensitivity to 95.0%.

8.1.1 Interpretation of results

Normal (negative)	< 2.3 U/mL
Cut-off	3.1 U/mL
Positive	> 3.9 U/mL
Equivocal	2.3 - 3.9 U/mL

The results alone should not be the only reason for any therapeutic consequences. The results should be correlated with other clinical observations and diagnostic tests.

9 LIMITATIONS OF USE

Reliable and reproducible results will be obtained, when the assay procedure is performed with a complete understanding of the instructions for use and with adherence to good laboratory practice. Any improper handling of samples or modification of this test might influence the results.

10 PERFORMANCE CHARACTERISTICS

10.1 Assay Dynamic Range

The dynamic range of the assay is defined by the limit of detection and the maximum value of the Master Curve.

Values found below the measuring range are indicated as " $< 0.179 \text{ U/mL}$ ".

Values found above the measuring range are indicated as " $> 100 \text{ U/mL}$ ".

The measuring range of the assay is between $0.179 \text{ U/mL} - 100 \text{ U/mL}$.

10.2 Sensitivity

The sensitivity study was designed according to CLSI guideline EP17-A2. The Limit of Blank (LoB) is 0.060 U/mL. The Limit of Detection (LoD) is 0.179 U/mL. The Limit of Quantification (LoQ) is 0.479 U/mL.

10.3 Precision Performance

The precision study was designed on the basis of to CLSI guideline EP5-A2.

10.3.1 Intra-Device

The Intra-Device precision was determined with 4 patient samples covering the measuring range in 5 independent runs on 5 days with 2 different devices in 5 replicates per run. CV was calculated as mean CV of 10 runs.

The within-run variability is shown below:

Sample	n	Mean (U/mL)	CV (%)
1	10	1.0	6.8
2	10	6.0	7.2
3	10	13.5	6.1
4	10	35.9	5.1

10.3.2 Inter-Device Precision

The Inter-Device precision was determined with 4 patient samples covering the measuring range in 5 independent runs on 5 days with 2 different devices with 5 determinations per run. CV was calculated from 50 determinations.

The total precision is shown below:

Sample	n	Mean (U/mL)	CV (%)
1	50	1.0	11.6
2	50	6.0	9.9
3	50	13.5	8.6
4	50	35.9	8.1

10.3.3 Inter-Lot

The between-lots variation was determined by 6 measurements of 4 samples with 3 different kit lots.

Sample	n	Mean (U/mL)	CV (%)
1	18	5.5	0.5
2	18	12.4	2.4
3	18	23.1	4.7
4	18	33.3	5.8

10.4 Recovery

Recovery was determined by adding four increasing concentrations of the analyte to four different patient samples containing different amounts of endogenous analyte. Each sample (non-spiked and spiked) was assayed and analyte concentrations of the samples were calculated from the Master Curve. The percentage recoveries were determined by comparing expected and measured values of the samples.

Sample	1	2	3	4
Concentration (U/mL)	6.4	13.8	33.7	51.6
Average Recovery (%)	100.2	96.5	92.4	103.9
Range of Recovery (%)	from 92.4	87.2	88.3	98.3
	to 105.2	100.2	98.4	110.4

10.5 Linearity

Four samples containing different amounts of analyte were serially diluted with Sample Diluent. The percentage recovery was calculated by comparing the expected and measured values for the analyte.

Sample	1	2	3	4
Concentration (U/mL)	33.5	47.9	86.1	97.7
Average Recovery (%)	104.3	94.0	106.3	107.8
Range of Recovery (%)	from 95.5	86.9	101.3	97.9
	to 108.6	100.7	114.9	114.5

10.6 Interfering Substances

Haemoglobin (up to 4 mg/mL), Bilirubin (up to 0.5 mg/mL) and Triglyceride (up to 7.5 mg/mL) have no influence on the assay results.


10.7 High-Dose-Hook Effect

Hook effect was not observed in this test up to a concentration of 2000 U/mL of anti-tTG IgA.

11 METHOD COMPARISON

A comparison of DRG:HYBRID-XL Anti-tTG IgA Test HYE-5894 (y) and the DRG t-Transglutaminase-Ab IgA ELISA (EIA-4259) (x) using clinical samples gave the following correlation:

$$\begin{aligned} n &= 115 \\ r &= 0.962 \\ y &= 0.76x + 0.865 \end{aligned}$$

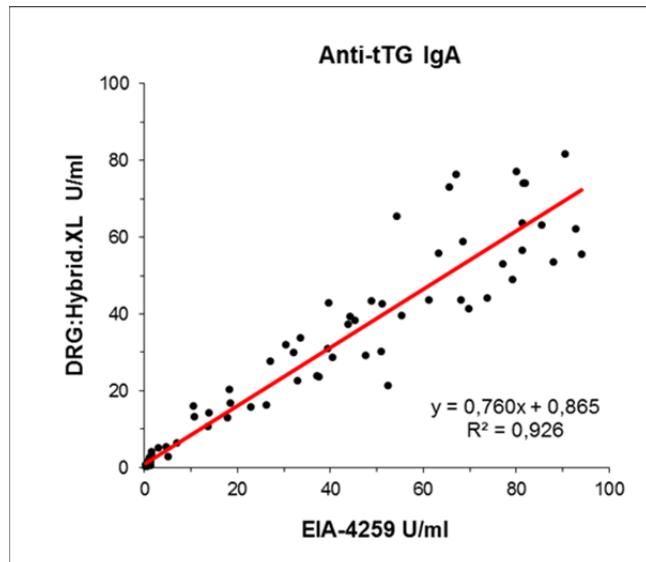12 LEGAL ASPECTS

12.1 Reliability of Results

The test must be performed exactly as per the manufacturer's instructions for use. Moreover, the user must strictly adhere to the rules of GLP (Good Laboratory Practice) or other applicable national standards and/or laws. This is especially relevant for the use of control reagents. It is important to always include, within the test procedure, a sufficient number of controls for validating the accuracy and precision of the test.

The test results are valid only if all controls are within the specified ranges and if all other test parameters are also within the given assay specifications. In case of any doubt or concern please contact DRG.

12.2 Therapeutic Consequences

Therapeutic consequences should never be based on laboratory results alone even if all test results are in agreement with the items as stated under point 12.1. Any laboratory result is only a part of the total clinical picture of a patient.

Only in cases where the laboratory results are in acceptable agreement with the overall clinical picture of the patient should therapeutic consequences be derived.

The test result itself should never be the sole determinant for deriving any therapeutic consequences.

12.3 Liability

Any modification of the test kit and/or exchange or mixture of any components of different kit lots could negatively affect the intended results and validity of the overall test. Such modification and/or exchanges invalidate any claim for replacement.

Claims submitted due to customer misinterpretation of laboratory results subject to point 12.2 are also invalid. Regardless, in the event of any claim, the manufacturer's liability is not to exceed the value of the test kit. Any damage caused to the test kit during transportation is not subject to the liability of the manufacturer.

For further information please refer to the User Manual of the DRG:HYBRID-XL, analyser-specific application sheets, product information and package inserts of all necessary components.

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1. Schuppan, D, & Zimmer, KP: The diagnosis and treatment of celiac disease. Deutsches Ärzteblatt International, 2013; 110(49): 835–46.
2. Sollid LM, Thorsby E: HLA susceptibility genes in celiac disease: genetic mapping and role in pathogenesis. Gastroenterology 1993; 105: 910-922.
3. Schuppan D, Junker Y, Barisani D: Celiac disease: from pathogenesis to novel therapies. Gastroenterology 2009; 137: 1912–33.
4. Fernández E et al.: Comparison of six human anti-transglutaminase ELISA-tests in the diagnosis of celiac disease in the Saharawi population. World J Gastroenterol 2005; 11(24):3762-3766.
5. Conrad, K, Schößler, W, Hiepe, F: Autoantikörper bei organspezifischen Autoimmunerkrankungen. Pabst Science Publishers, 2010; 309-314.
6. Herold, G. Innere Medizin. 2016; 477-478.
7. Lorand, L, Conrad, SM:Transglutaminases. Mol Cell Biochem, 1984; 9-35.
8. Di Sabatino, A et al.: The function of tissue transglutaminase in celiac disease. Autoimmunity Reviews; 2012.
9. Dieterich, W et al.: Identification of tissue transglutaminase as the autoantigen of celiac disease. Nature Medicine, 1997; 3(7): 797–801.
10. Krause, I et al.: Anti-infectious antibodies and autoimmune-associated autoantibodies in patients with type I diabetes mellitus and their close family members. Annals of the New York Academy of Sciences. 2009; 1173: 633–9.
11. Farrace, MG et al.: Presence of anti-tissue transglutaminase antibodies in inflammatory intestinal diseases: an apoptosis-associated event? Cell Death & Differentiation.2001; 8 (7): 767–70.
12. Teichmann, J, Voglau MJ, Lange U: Antibodies to human tissue transglutaminase and alterations of vitamin D metabolism in ankylosing spondylitis and psoriatic arthritis. Rheumatology International. 2009; 30 (12): 1559–1563.
13. Picarelli, A et al.: Anti-tissue transglutaminase antibodies in arthritic patients: a disease-specific finding? Clinical Chemistry, 2003; 49 (12): 2091–4.
14. Fleckenstein B et al.: Molecular characterization of covalent complexes between tissue transglutaminase and gliadin peptides. J. Biol. Chem. 2004; 279 (17): 17607–16.
15. Thomas, L.: Antikörper bei gastrointestinalen Erkrankungen. In: Thomas, L. (ed.): Labor und Diagnose (8th edition, 2012), 1508 - 1514, TH-Books Verlagsgesellschaft mbH, Frankfurt/Main
16. Mäki M, et al.: Prevalence of Celiac disease among children in Finland. N Engl J Med 2003; 348: 2517-2524.
17. Ventura A, Magazza G, Greco L: Duration of exposure to gluten and risk for autoimmune disorders in patients with celiac disease. SIGEP Study Group for Autoimmune Disorders in Celiac Disease. Gastroenterology 1999; 117: 297–303.18.

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